

Calculator Smoke 1

1. Concentration

- Millipore: Human, Mouse, Rat: 20X for Beads & ABS. 1X on PE.
- Thermofisher: Human, Mouse, Rat: 50X for Beads & ABS. 1X on PE.
- Bio-Rad: Human, Mouse, Rat: Variable, 20X or 10X (analyte dependent). 100X on PE.

2. Volume

- Millipore: Bead, Antibodies (ABS), and PE: 25 uls per well.
- Thermofisher: Beads & PE: 50 uls per well, Antibodies: 25 uls per well.
- BioRad: Beads & PE; 50 uls per well, Antibodies: 25 uls per well.

3. Input

- Concentration must be on a table with Analyte.
- Volume of Bead, ABS, PE as a variable input depending on Platform, Species, and Panel.
- Concentration and volume must be variable and can be easily adjusted by manager.
- These values should be in a table or datasheet most likely to excel.

4. Selection parameter

- Selections should be limited to per Panel.
 - Example: If Millipore Mouse Panel I is selected, therefore only analytes in Millipore Mouse Panel I should be available for selection.
 - Panel I and Panel II should not be an option for selection together.

Select Analytes

Premix Panels

No Premix (Custom Assay)
Select individual analytes only

Selected

Milliplex Mouse Premix Panel I 32-Plex
Comprehensive 32-plex cytokine panel for mouse samples

Milliplex Mouse Premix Panel II 23-Plex
23-plex cytokine panel for mouse samples

MILLIPLEX MOUSE PREMIX PANEL I 32-PLEX

IL-1a
32 20x

IL-1b
47 20x

IL-2
58 15x

IL-6
19 25x

IL-17
47 20x

GM-CSF
22 15x

TNF-a
75 20x

MILLIPLEX MOUSE PREMIX PANEL II 23-PLEX

IFN-g
33 20x

IL-4
53 20x

IL-10
18 20x

Selected Analytes (3)


IL-1a 32

GM-CSF 22

IL-4 53

5. Duplicate error not in proper sequence. It should be full column not row.

Plate Layout

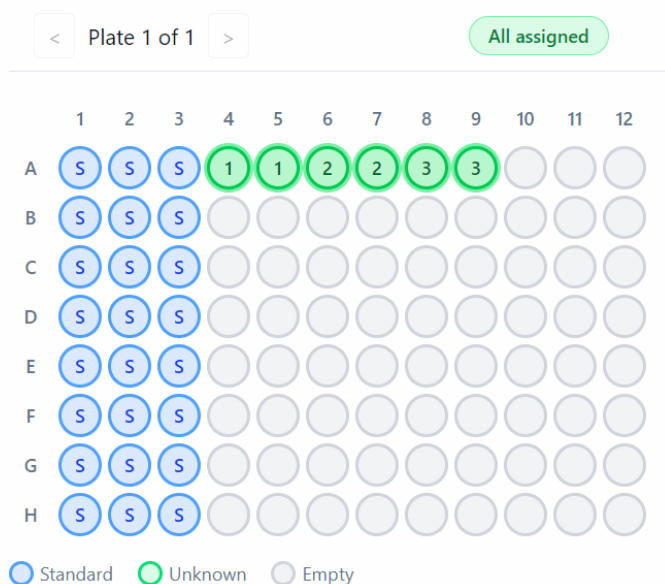

6. Error on adding new plate.
- Adding new plate removes previous plate record.
7. Summary input
- Species did not carry over in the summary input
 - Panel still should show what panel the analyte selected from.
8. Summary Out
- Include Summary Input
 - Include Actually calculation work
 - PDF or printable version